## LifeKit by MTTS



# LIGHTMETER

An accurate instrument to determine the intensity of phototherapy light.

#### Overview

- Model LM-800 - compact, handheld unit

#### Sensor

- spectral response chosen to match wavelength used for
- jaundice treatment
- two sided measurement for convenient measurement
- of double sided phototherapy devices such as the MTTS Firef y


#### Displays/Indicators

- touch screen liquid-crystal display (LCD)
- alarms: low battery, system failure, out of measurement
- results 'freeze'
- bandwidth estimator
- absolute and relative irradiance measurement
- multi-point mode

#### Additional Features - available in 2022

- oxygen measurement
- f ow measurement

### Spectral Response

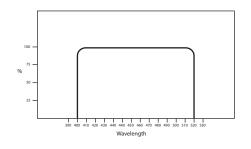



#### Dimensions and Weight

Height 148 mm Width 73 mm Depth 17 mm Weight 200 g

#### Performance Specif cations

Application LED and blue f uorescent phototherapy light sources; using the lightmeter for halogen phototherapy results

in a minor but not considerable error

400-520nm Spectral response

Measurement range single side

 $0.1 - 150.0 \,\mu W/cm^2/nm$  $0.1 - 300.0 \,\mu W/cm^2/nm$ double side  $0.1 \,\mu W/cm^2/nm$ Resolution

±2% @ 30 degree angle Cosine characteristics ±7% @ 60 degree angle

±25% @ 80 degree angle +/- 3% of reading

#### **Battery Specif cation**

Accuracy

Type 2xAAA 1.5V Alkaline, Non-rechargeable Lifetime Up to 120h operating time

#### **Environmental Specif cation**

Operation Temperature  $10^{\circ}\text{C}$  to  $40^{\circ}\text{C}$ Storage Temperature  $0^{\circ}$ C to  $50^{\circ}$ C Operation Atmospheric Pressure 70 - 106 kPA Storage Atmospheric Pressure 70 - 106 kPA

Operation Humidity 30 – 90% RH non condensing Storage Humidity 5 - 90% RH non condensing

Sealing

Drop/Shock 1m onto hard surface, meets or exceeds IEC 60068-2-32

(only with the protective case)

#### Standards

ISO 13485:2016 IEC 61010-1:2010 IEC 61326-1:2012

#### Ordering Information and Servicing

ORDERING INFO ITEM ORDER NO Lightmeter LM-800

Rechargeable batteries LM-800-010

and charger

Calibration

Parts

Required

Every 2 years after delivery date; return to MTTS for calibration

No customer serviceable parts; return to MTTS for service

Main Unit

Optional

MEDICAL TECHNOLOGY TRANSFER AND SERVICES Co., LTD ISO 13485: 2016 Certif ed Company

#### Distributeur:



@ info@ckapaca.com

www.ckapaca.com

**6** +39 349.805.29.90

32100 BELLUNO, Italie